
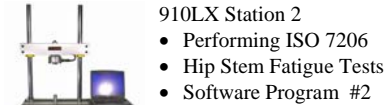


2370 MS Laboratory Management System for Medical Device Test Laboratory – 6 Channel Example

- 910LX Servohydraulic Fatigue Test System – 2 stations shown
- 831L Axial Torsion Test System (2 channels) – 1 station
- 800LM Low Force High Frequency Fatigue – 1 station
- 310 Series Electromechanical (Static) Test Machine – 1 station

Each Test System includes

- Digital Display, Keypad, Mouse (not shown)
- Network Connection to 2370MS – No PC
- 24 Bit Signal Conditioners:
 - Load
 - Position
 - Strain
 - Channel 4 (open)

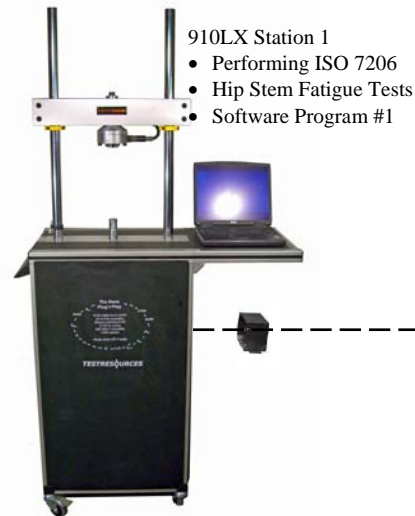

2370MS

- Max 16 servo-output channels
- Max 80 feedback channels
- 48 Digital output bits
- 32 Digital input bits



MS Control Software

- Central Lab Data Storage
- Monitors Software #1-#5
- Global Data Sharing (GDS)

